FarmaMed

**MASSAGGIATORE ADDOME E FASCIA LOMBARE
LUMBAR AND ABDOMINAL MASSAGER
MASAJEADOR LUMBAR Y ABDOMINAL
MASSEUR ABDOMINAL ET LOMBAIRE
BAUCH- UND LENDEN MASSAGEGERÄT
MESSAGEADOR ABDOMINAL E CINTA LOMBAR**

ISTRUZIONI PER L'USO	IT
INSTRUCTIONS FOR USE	EN
INSTRUCCIONES DE USO	ES
INSTRUCTIONS D'UTILISATION	FR
GEBRAUCHSANLEITUNG	DE
INSTRUÇÕES DE USO	PT

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LUMBAR AND ABDOMINAL MASSAGER

Dear customer, Gabbiano S.p.A. wishes to thank you
for choosing this **Farmamed** product.

Our brand is synonymous with reliable, high-quality products.

For further information on this or other **Farmamed** products,
visit our website [**www.farmamed.it/en/**](http://www.farmamed.it/en/)

INDEX

PAG

- 20 _ Package contents
- 20 _ Description of the appliance
- 22 _ Symbols key
- 22 _ Security warning
- 25 _ Warnings on rechargeable battery
- 25 _ Warnings on disposable battery
- 26 _ Warnings on electromagnetic compatibility
- 27 _ Appliance preparation
- 28 _ Instruction for use
- 30 _ Remote control battery replacement
- 30 _ Maintenance and upkeep
- 30 _ Disposal
- 31 _ After-sales service and warranty
- 32 _ Trouble-shooting
- 33 _ Product specifications

PACKAGE CONTENTS

- Massager
- Remote control with inserted CR2032 battery
- USB cable
- This user manual

Device functions

This device generates electrical impulses that help contract and passively relax the muscle fibers of the abdominal and lumbar wall. This stimulation helps warm up and relax the muscles, as well as define and tone them (although it cannot be considered a substitute for active training). It also alleviates the feeling of fatigue and promotes blood circulation. Thanks to the integrated rechargeable battery, the device can be used anytime and anywhere.

Intended for personal and domestic use only.

It is not a medical device.

Device features

- Ergonomic design, in line with the curve of the waist.
- Remote control of the device, thanks to the included remote control with an LCD display.
- Simple and intuitive menu on the remote control for navigation.
- Various functions in one device: infrared hot plate, low-frequency electrical pulse, combinable through 6 programs.
- Possibility to select the desired temperature level and the intensity of the electrical pulse according to your needs and preferences.

DESCRIPTION OF THE APPLIANCE

Fig. 1: main unit

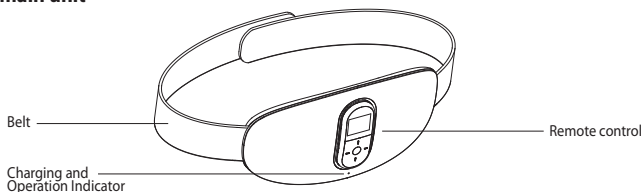
Fig. 2: main unit, interior view

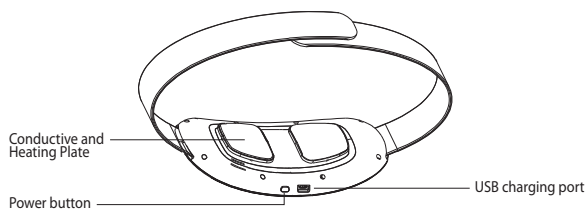
Fig. 3: remote control, front


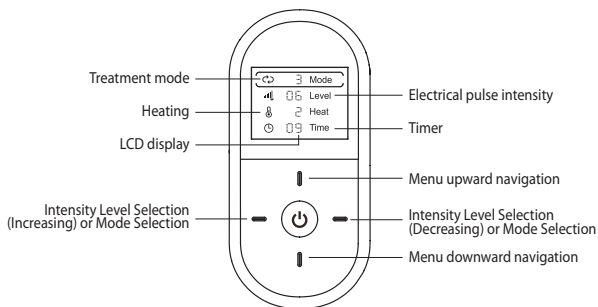
Fig. 4: remote control, back

SYMBOLS KEY

The following symbols appear in the instructions for use and on the appliance and its packaging.

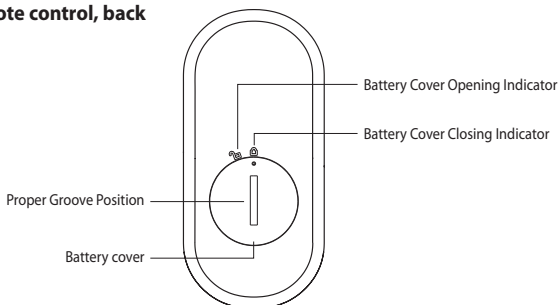Warning – Warning regarding possible risk of injury or possible health hazards.

Note – Note providing important information.

Dispose of exhausted and completely discharged batteries in the appropriate collection points, at waste collection site for toxic waste or at electronics stores. Battery disposal is a legal requirement.

Disposal in compliance with the provisions set forth in Directive 2012/19/EU on waste electrical and electronic equipment (WEEE).

Serial number.

Please read this manual carefully before using.

According to the degree of protection against direct and indirect contacts, the following equipment is classified as a bf type.

The CE mark certifies compliance with the safety and health requirements outlined in the relevant directives for electrical and electronic products and components, falling under the Low Voltage Directive 2014/35/EU, Electromagnetic Compatibility Directive (2014/30/EU), and RoHS (2011/65/CE) + AEE (2015/863/UE).

IP20 Degree of protection of electrical equipment covers, where the first number indicates the degree of protection against penetration by solid foreign bodies (from 0 to 6) and the second number the degree of protection against liquid penetration (from 0 to 8).

SECURITY WARNING

Before using the appliance, read the instructions for use carefully, paying particular attention to the warnings regarding safety, which must be observed. This manual must be considered part of the product and must be kept for its entire life-cycle. If the appliance is passed on to another owner, it must be provided together with all the relevant documentation. If the manual is mislaid or for further information or clarification, please contact Gabbiano S.p.A. using the details provided on the last page.

- Before using the product, check that it is intact and has no visible damage. If you have any doubts, do not use the appliance and contact your retailer.

- Using the electrode close to the chest may increase the risk of atrial fibrillation.
- Do not dismantle, repair or tamper with the device. Avoid it.
- If you feel unwell or experience abnormal physical and/or physiological reactions during use, discontinue use of the product and seek medical advice immediately.
- Do not use the appliance on areas close to the heart, head, eyes, front of the neck (in particular, the carotid artery), the lower back, oral cavity or genitalia.
- Do not use on areas affected by skin conditions.
- The device must be switched off and switched back on again before using it, if it is moved or when changing the area being treated.
- If you feel unwell while using the product, discontinue use immediately and seek medical advice.
- Do not use together with other electronic medical devices such as pacemakers, artificial heart or lungs and other electronic medical devices with lifesaving systems, ECG machines and other electronic medical devices.
- Do not use the product in places subject to high temperatures, in the presence of flammable substances or in areas subject to electromagnetic radiation and humidity.
- In order to avoid risks, the main unit's battery must be replaced by a professional maintenance technician.
- The skin should not be treated with oily creams or ointments before use.
- Keep the appliance away from heat sources and do not use it in the vicinity (~1 m) of short-wave or microwave appliances (e.g. mobile telephones), in order to avoid undesirable current surges.
- Do not expose the appliance to direct sunlight or high temperatures.
- Protect the appliance from dust, dirt and humidity. Do not immerse the appliance in water or other liquids.
- This appliance is intended for personal use, not for a professional one.
- Remove any metal objects (such as belt buckles or necklaces) that could come into contact with the electrodes during use. Any jewellery or piercings present in the area to be treated (e.g. navel piercing) should be removed before using the appliance.
- Keep the appliance out of the reach of children, do not use it on them unless under constant medical supervision.
- When the device is moved, it must be handled with care, avoiding sudden shaking.
- Pay attention to the device protection, avoiding abrasion.

Limitations

The appliance is intended solely for the treatment therapy of the abdominal and lumbar area of the human body. It must not be considered a substitute for medical treatment. If you are unsure whether use of the appliance is appropriate, contact your doctor.

Use of the device is forbidden in the following cases:

- patients with pacemakers or other electrical implanted devices and those with suspected or confirmed epilepsy;
- in the presence of wet skin or perspiration or when sleeping;
- patients with brain haemorrhages;
- pregnant women;
- in the presence of a haemorrhage following serious injury or a fracture;
- subjects with purulent inflammation, blood poisoning;
- in the presence of medical conditions or lesions in the area of the body to be treated;
- subjects with acute cardiovascular or cerebrovascular disease;
- individuals with sensitive skin, heart disease, abnormal blood pressure, malignant tumours, cerebrovascular disease, acute medical conditions or who are undergoing therapy or have other medical conditions;
- children or individuals unable to use the device themselves.

Consult your doctor before using the device in the presence of:

- a serious medical condition or after upper-body surgery or subjects who suffer from thrombosis or diabetes;
- chronic pain of unknown origin, regardless of the area of the body;
- acute medical conditions, especially in the case of suspected or confirmed diseases associated with hypertension, blood clotting disorders, history of thromboembolic disease or cancer;
- any kind of skin disease or open wounds;
- swelling, skin infections, skin disease or phlebitis;
- skin sensitivity disorders or in those who are not sensitive to heat or electrical stimulation applied to the skin;
- on-going medical treatment;
- disorders that present during stimulation treatment;
- skin irritation caused by prolonged stimulation in the same point.

The product is made of non-surgical quality stainless steel in parts intended to come into direct contact with the skin. In case of doubts about the possible onset of sensitization or allergies, consult your doctor before using the device. **Stainless steel plates contain a small amount of nickel.**

WARNINGS ON RECHARGEABLE BATTERY

The device has a built-in battery which is not accessible to the user and which can be recharged using the USB cable supplied (power supply not included). It cannot be replaced by the user, who must contact Gabbiano S.p.A.

- Turn off the device when connecting the USB cable for charging.
- Do not turn on the device when the USB cable is inserted.
- Do not use the unit while recharging.
- Use only the micro-USB cable provided by the manufacturer.
- Use only the chargers for the main unit specified in the "Product specifications" section of this manual.
- Before first use, the unit must be charged. When using it for the first time, fully charge the rechargeable battery for at least 30 minutes.
- **Danger of ingestion!** If the device is broken, children could swallow the rechargeable battery and suffocate.
- **Risk of explosion!** DO NOT throw the device with the battery into a fire, even if it is exhausted.
- Do not disassemble, open or crush the internal battery.
- If the battery of the main unit leaks, wear protective gloves and dispose of the device.
- If the liquid from the rechargeable battery comes into contact with your skin or eyes, rinse the affected areas with water and seek medical attention.
- The rechargeable battery must be correctly charged before use. Always follow the manufacturer's instructions and the information provided in this user manual to charge the battery correctly (see section "Preparing the device").
- To achieve the longest possible life cycle, fully charge the battery pack at least once a month.

WARNINGS ON DISPOSABLE BATTERY

If the appliance is not to be used for a long period, remove the batteries and store in a cool, dry place at room temperature.

- DO NOT recharge batteries unless they are rechargeable.
- Rechargeable batteries must ONLY be recharged as indicated in the manual provided with them and used for the appliances indicated by the manufacturer.

- NEVER expose batteries to sources of heat and direct sunlight. Failure to observe this precaution may damage the batteries and/or cause them to explode.
- Do NOT incinerate, even when the batteries are exhausted. Risk of explosion!
- Use batteries of the same type and replace them at the same time.
- Batteries must be removed or replaced by an adult.
- Keep the batteries out of the reach of children: swallowing batteries can be fatal. If swallowed, seek medical advice immediately.
- The acid inside the batteries is corrosive. Avoid contact with skin, eyes and clothing. If the battery liquid comes into contact with the skin and eyes, wash the affected area with water and seek medical advice.
- In the case of battery liquid leakage, use protective gloves and clean the battery compartment with a dry cloth. Dispose of the materials used for cleaning.

WARNINGS ON ELECTROMAGNETIC COMPATIBILITY

The appliance is suitable for use in any condition indicated in these instructions for use, including domestic environments.

- Under certain circumstances, in the presence of electromagnetic disturbance, use of the appliance may be subject to limitations. Consequently, error messages may appear on the display/appliance/remote control.
- Do not use this appliance in the immediate vicinity of other electrical appliances, as this could cause a malfunction.
- Failure to comply with these notes could be detrimental to the performance of the appliance.
- Mobile telephones or appliances that generate electromagnetic fields must not be used in the vicinity of the product.
- Do NOT use the product in environments subject to electromagnetic interference.
- Do NOT use simultaneously with other devices that transmit electrical impulses.

EN

26

Additional quantitative data regarding electromagnetic compatibility are available from Gabbiano S.p.A. and can be obtained upon request.

APPLIANCE PREPARATION

INSTALLATION OF THE REMOTE CONTROL BATTERY


When the brightness of the LCD display on the remote control diminishes or the device does not respond to commands, it is likely that the remote control battery is running low, and you need to replace it with a new one.

The remote control uses a CR2032 battery (included in the package).

To insert the battery, open the battery compartment (Fig. 5).

Use a coin or a flathead screwdriver in the groove to unscrew the battery cover, rotate it to "⌚" to loosen the cover. Insert the battery while observing the indicated polarity. Close the cover by aligning the battery cover mark with "⌚" and screw it towards "🔒".

Fig. 5

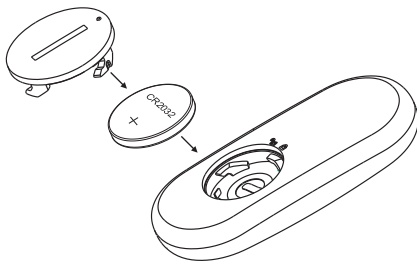INTEGRATED RECHARGEABLE BATTERY

The device is equipped with a built-in rechargeable battery that is not user-accessible, allowing for usage and operation without external power.

- 1) During device use, if the lithium battery charge level is insufficient, the device emits an audible signal 3 times. To recharge the battery, insert the USB C connector of the included cable into the USB charging port of the device, then connect the USB cable connector to an appropriate power source (not included, DC 5V) compliant with electrical safety standard IEC 60601-1. At this point, the charge indicator flashes. The charging time is approximately 2/3 hours. Once the battery is fully charged, the indicator remains on.
- 2) When rechargeable batteries are not used for an extended period, their lifespan decreases. It is recommended to perform a full recharge at least once a month.

INSTRUCTION FOR USE

Introduction

It is recommended to use the device a maximum of 2 times a day. During each session, use the electrostimulation function for no more than 10-15 minutes and the heating function for no more than 15-30 minutes.

- 1) Wear the massager in a way that the metal plates adhere well to the skin (do not use it over clothing). Adjust the belt so that it is neither too tight nor too loose.
- 2) Press the device's power button: you will hear an audible signal, and the operation indicator will start flashing.
- 3) Press the "U" button to turn on the remote control: you will hear an audible signal, the screen backlight will illuminate, the operation indicator will stop flashing, and the device will start functioning. The following default settings will be displayed on the screen:
 - Treatment mode: 1;
 - Electrical pulse intensity: 00;
 - Medium heating: 2;
 - Timer: 15 minutes.

To stop the operation, press the power on/off button on the remote control at any time.

Note: If no operation is performed within 10 seconds after turning on, the remote control will enter power-saving mode, and the screen brightness will decrease to conserve battery capacity. Pressing the remote control will restore the normal screen brightness.

Explanation of the Wireless Remote Control and Device Functions

The remote control features the following buttons:

"↑" Button: Press the top button to select a menu item by moving upward and the bottom button to select a menu item by moving downward.

"←" Button: Press the left button to decrease the intensity level or the selected mode, and press the right button to increase the intensity level or the selected mode.

- **Treatment Mode “MODE”** ().

Select “” using the “” buttons, then press “” to choose one of the 4 modes of the device:

- 1) **Automatic Mode 1:**

Combines massage modes 1-2-3, repeating them in cycles

- 2) **Automatic Mode 2:**

Combines massage modes 4-3-6, repeating them in cycles

- 3) **Massage Type 1**

- 4) **Massage Type 2**

- 5) **Massage Type 3**

- 6) **Massage Type 4**

After selecting the mode, you can adjust the remaining parameters using the provided wireless remote control:

- **Electrical Pulse Intensity “LEVEL”** ().

Select “” using the “” buttons, then press “” to increase/decrease the electrical pulse intensity. There are 16 ascending levels of electrostimulation intensity.

- **Heating Level “HEAT”** ().

Select “” using the “” buttons, then press “” to choose the desired heating level from:

“0”: No heating

“1”: Low heating

“2”: Moderate heating

“3”: High heating

The default setting at startup is moderate heating.

- **Adjustment of Treatment Duration “TIME”** ().

Select “” using the “” buttons, then press “” to select the treatment duration. The adjustable time interval ranges from 5 to 30 minutes (in 5-minute intervals).

REMOTE CONTROL BATTERY REPLACEMENT

- 1) Open the battery compartment as indicated in figure 5.
- 2) Remove the exhausted batteries and dispose of them as indicated in the "Disposal" chapter.
- 3) Insert batteries of the same type and with the same service life, making sure they are positioned as indicated.
- 4) Dispose of the exhausted batteries in compliance with applicable local regulations.
- 5) Close the battery compartment.
 - The remote control is powered by a CR2032 battery.
 - If the product is not to be used for a long period, remove the battery, which could otherwise leak and damage the device.
 - Do not reverse the positive and negative poles of the battery; the positive pole must face upwards.
 - Dispose of exhausted batteries in compliance with applicable local regulations. See the "Disposal" chapter of this manual.

MAINTENANCE , UPKEEP AND CLEANING

After each use, wipe the conductive plate with a damp soft cloth and dry.

- If the conductive plate is particularly dirty, clean with a soft cloth wetted with a small amount of rubbing alcohol (75% solution).
- The device must be kept in a dry, well-aired and safe place.
- The device must be transported with care, avoiding excessive shaking.
- It should be adequately protected to avoid scratching.

DISPOSAL

Dispose of exhausted and completely discharged batteries in the appropriate collection points, at waste collection site for toxic waste or at electronics stores. Battery disposal is a legal requirement.

For environmental reasons, the appliance must not be discarded with normal household waste at the end of its service life. It must be taken to an appropriate disposal facility. Dispose of the appliance in accordance with European Directive 2012/19/EU regarding waste electrical and electronic equipment. For further information or clarification, contact the local authority responsible for waste disposal.

The European Directive 2012/19/EU on Waste Electrical and Electronic Equipment stipulates that these devices should not be disposed of in the normal flow of urban solid waste but should be collected separately to optimize the recovery and recycling of the materials they contain and to prevent potential health and environmental hazards due to the presence of potentially dangerous substances. Waste can be taken to designated collection centers or can be returned free of charge to the retailer when purchasing a new equivalent device, or without the obligation to make a purchase for devices smaller than 25 cm. Sanctions are provided for illegal disposal. For further information on the proper disposal of these devices, you can contact the relevant public service.

AFTER-SALES SERVICE AND WARRANTY

This appliance is intended for domestic use only and is not suitable for professional purposes.

- The appliance is covered by a two-year warranty from the time of purchase, which must be proven by proof of purchase. This period complies with applicable regulations ("Consumer Code" Italian Legislative Decree no. 206 of 06/09/2005) and only applies if the consumer is a private subject.
- The warranty covers production defects only and is not valid for damage caused by accidents, incorrect use, negligence or improper use of the product. It does not cover materials subject to normal wear and tear.
- The appliance must not be opened under any circumstance; the warranty is not valid if the product has been opened or tampered with.
- The warranty expires 2 years after purchase; after this time, any technical assistance will be provided on a paid basis. For information on technical assistance under warranty or provided on a paid basis, contact Gabbiano S.p.A. using the details on the last page of this manual. Consumers will not receive any reimbursement for product

repairs or replacements performed under warranty. In the event of a malfunction, contact the manufacturer; Do NOT send anything directly to Gabbiano S.p.A., **send an e-mail to service@farmamed.it first**. Work performed under warranty (including replacement of all or part of the product) will not extend the duration of the original warranty of the replaced product.

- The manufacturer declines all responsibility for any damage/injury that may be directly or indirectly caused to people, property or domestic animals as a consequence of failure to comply with the indications provided in this instruction manual, in particular the warnings regarding the installation, use and maintenance of the appliance.
- **Given its on-going commitment to improving its products, Gabbiano S.p.A. reserves the right to modify all or parts of its products in order to meet production requirements, without giving advanced notice, without this entailing any responsibility for Gabbiano S.p.A. or its dealers.**

TROUBLE-SHOOTING

Please check the following points if you encounter any issues during the use of the device.

Problem	Possible cause	Solution
The remote control doesn't work.	The battery is dead.	Replace it with a new CR2032 battery.
The remote control works, but the device doesn't turn on.	The internal battery's charge level is insufficient.	Recharge the battery using an external power source.
	The device is connected to the power supply, but the internal battery level is insufficient.	Ensure that the device is properly connected to the power supply, and verify that the power supply is functioning correctly.
Unable to recharge or no charging indication	The adapter is not connected correctly.	Check if the power adapter is properly connected to the device and ensure there are no abnormal conditions (e.g., wire breaks or bends).
	The USB cable is not properly to the device.	Verify if both ends of the USB cable are correctly inserted.
The device turns on but doesn't seem to work	The electrodes are not in contact with the skin.	Ensure that the electrodes are in direct contact with the skin.
	The pulse intensity has not been adjusted.	Once the device is switched on, the initial intensity is set to 0, so it needs to be adjusted according to your needs.

PRODUCT SPECIFICATIONS

Product name: Lumbar and abdominal massager

Product dimensions (without belt): 26 x 13 x 2.5 cm

Weight: approximately 390 g

Product power supply requirements

- **Rated power:** 5 W
- **Internal electrical power supply:** 3.7 V DC. 2000mAh
- **Remote control power supply:** CR2032 3 V DC

Environmental product requirements

Operating conditions:

- **Temperature:** 5 °C - 40.0 °C
- **Maximum relative humidity:** ≤ 80% RH
- **Atmospheric pressure:** 86 - 106 kPa

Storage and transportation conditions:

- **Temperature:** -20 °C - +55 °C
- **Maximum relative humidity:** ≤ 93% RH
- **Atmospheric pressure:** 50.0 - 106.0 kPa

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