# FarmaMed

**TERMOMETRO A INFRAROSSI SENZA CONTATTO**  
**INFRARED NON-CONTACT THERMOMETER**  
**TERMÓMETRO INFRARROJO SIN CONTACTO**  
**THERMOMÈTRE INFRAROUGE SANS CONTACT**  
**KONTAKTLOSES INFRAROT THERMOMETER**  
**CONTACTLOZE INFRAROOD THERMOMETER**  
**BERÖRINGSFRI INFRARÖD TERMOMETER**  
**TERMÓMETRO DE INFRAVERMELHOS SEM CONTACTO**

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**Schermo:** schermo LCD

## INFRARED NON-CONTACT THERMOMETER

Dear customer, **Farmamed** wishes to thank you for choosing this product.

Our name is synonymous with reliable, high-quality products.

For further information on this or other **Farmamed** products, visit our website  
**[www.farmamed.it](http://www.farmamed.it)**.

### INDEX

#### PAG

- 19 \_ Package contents
- 19 \_ Description of the appliance
- 20 \_ Symbol key
- 21 \_ Security warning
- 22 \_ Warnings on battery use
- 23 \_ Warnings on electromagnetic compatibility
- 24 \_ Appliance preparation
- 25 \_ Basic parameters
- 26 \_ Instructions for use
- 28 \_ Battery replacement
- 29 \_ Maintenance and upkeep
- 29 \_ Disposal
- 30 \_ After-sales service and warranty
- 31 \_ Trouble-shooting
- 32 \_ Product specifications

## PACKAGE CONTENTS

- Thermometer (replaceable batteries provided)
- These instructions for use

### Appliance Functions

This appliance can be used to measure human body temperature on the forehead without skin contact.


It can also be used to measure the temperature of liquids/objects and environments.

### Characteristics of this appliance


- Back-lit LCD display in three colours (green, orange and red), for an immediate interpretation of the body temperature measurement.
- One-second temperature detection.
- Memorises up to 9 measurement.
- Switches off automatically after 30 seconds' inactivity.
- Silent temperature detection.

## DESCRIPTION OF THE APPLIANCE

**Fig. 1: view of the appliance**

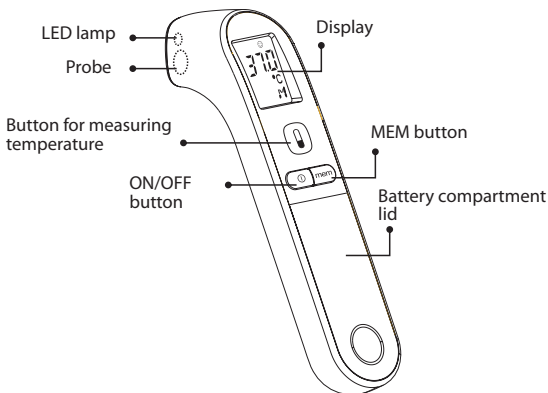**Fig. 2: display**
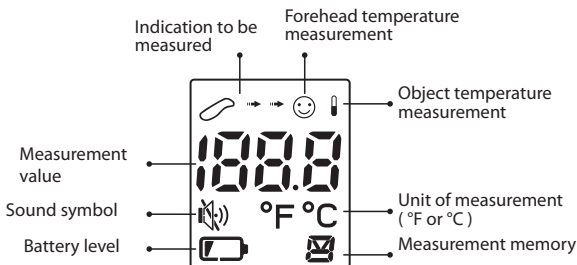
## SYMBOL KEY

The following symbols appear in the instructions for use and on the appliance and its packaging.

**Warning** – regarding the risk of injury or health hazards.

**Nota** – Note providing important information.

BF type applied part.

Protect from humidity.

Read the manual carefully before use.

Dispose of in compliance with the provisions of Directive 2012/19/EU regarding waste electrical and electronic equipment.

Dispose of exhausted and completely discharged batteries in the appropriate collection points, at waste collection site for toxic waste or at electronics stores. Battery disposal is a legal requirement.

CE label certifies that the device is compliant with the basic requirements set out in the legislation concerning the safety on medical devices.

**IP22:** Degree of protection of electrical equipment covers, where the first number indicates the degree of protection against penetration by solid foreign bodies (from 0 to 6) and the second number the degree of protection against liquid penetration (from 0 to 8).

Serial number.

## SECURITY WARNING

Before using the appliance, read the instructions for use carefully, paying particular attention to the warnings regarding safety, which must be observed. This manual must be considered part of the product and must be kept for its entire life-cycle. If the appliance is passed on to another owner, it must be provided together with all the relevant documentation. If the manual is mislaid or for further information or clarification, please contact Farmamed using the details provided on the last page.

- Before using the product, check that it is intact and has no visible damage. If you have any doubts, do not use the appliance and contact your dealer.
- This appliance must be utilised for its intended use only and as described in the instructions for use. Any other use is to be considered improper and, therefore, dangerous. The manufacturer declines all responsibility for any damage caused by incorrect or improper use.
- The measurements obtained using this unit must be considered as a guide and may not, under any circumstance, be considered as a substitute for diagnosis by a medical practitioner. Consult your doctor for explanations on the body temperature measured and evaluate any therapies on his indication.
- Keep out of the reach of children. Children of 12 years or under must use this unit under adult supervision.
- Do not touch or blow on the infrared sensor to avoid distorting the measurement. If the infrared sensor is damaged or soiled, the measurements may be unreliable.
- In the event of a fault and/or malfunction, switch the appliance off without tampering with it. Always contact your dealer when repair work is required.
- Do NOT try to dismantle, repair or replace parts of this device yourself. Always contact your dealer.
- Treat the appliance with care, protect it from blows, extreme temperature changes, humidity, dust, direct sunlight and sources of heat.
- Do NOT allow liquids to penetrate inside the body of the thermometer. Never immerse the thermometer in water. Never rinse under running water. This appliance is not waterproof.

## WARNINGS ON BATTERY USE

- If the appliance is not to be used for a long period, remove the batteries and store in a cool, dry place at room temperature.
- DO NOT recharge batteries unless they are rechargeable.
- Rechargeable batteries must ONLY be recharged as indicated in the manual provided with them and using the appliances indicated by the manufacturer.
- NEVER expose batteries to sources of heat and direct sunlight.
- Failure to observe this precaution may damage the batteries and/or cause them to explode.
- Do NOT incinerate, even when the batteries are exhausted.

### **Risk of explosion!**

- Use batteries of the same type and replace them at the same time.
- Batteries must be removed or replaced by an adult.
- Keep the batteries out of the reach of children: swallowing batteries can be fatal. If swallowed, seek medical advice immediately.
- The acid inside the batteries is corrosive. Avoid contact with skin, eyes and clothing.
- If the battery liquid comes into contact with the skin and eyes, wash the affected area with water and seek medical advice.
- In the case of battery liquid leakage, use protective gloves and clean the battery compartment with a dry cloth. Dispose of the materials used for cleaning.

## WARNINGS ON ELECTROMAGNETIC COMPATIBILITY

- The appliance is suitable for use in any environment indicated in these instructions for use, including domestic environments.
- Under certain circumstances, in the presence of electromagnetic disturbance, use of the appliance may be subject to limitations. Consequently, error or failure messages may appear on the display or appliance.
- Do not use this appliance in the immediate vicinity of other electrical appliances, as this could cause a malfunction.
- Failure to comply with these notes could be detrimental to the performance of the appliance.
- Mobile telephones or appliances that generate electromagnetic fields must not be used in the vicinity of the product.
- Do NOT use the product in environments subject to electromagnetic interference.

Additional quantitative data regarding electromagnetic compatibility are available from Farmamed and can be obtained upon request.


## APPLIANCE PREPARATION

### 1) INSERTING THE BATTERIES


This appliance operates on two replaceable AAA 1.5 V alkaline batteries (provided).

- To insert the batteries, open the battery compartment by pressing down on the lid, which will be released automatically (Fig. 3).
- Insert the batteries as indicated, then close the battery compartment.
- When the battery is exhausted, a message will appear on the display (Fig. 4).

**Fig. 3**

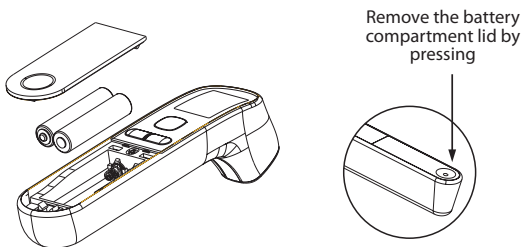**Fig. 4**
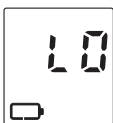
### 2) SELECT MEASUREMENT UNIT

- The default setting of this appliance is Celsius degrees.
- The scale of measurement can be switched between Fahrenheit degrees (°F) and Celsius degrees (°C).
- To change the measurement unit, while the appliance is switched off, hold down the "MEM" button.
- Change the measurement unit by pressing the "MEM" button until the display shows the desired unit of measurement.
- Wait for a few seconds for the measurement unit to be memorised. The appliance will switch on automatically.

## BASIC PARAMETERS


### 1) SOUND FUNCTION

- When the appliance is switched on, press the “MEM” key to enable or disable the sound.
- Press the “MEM” button, the “🔊” symbol will appear on the LCD display; this means that the sound is on. In the meantime, you will hear a short tone.
- Press the “MEM” key again: the “🔊” symbol will disappear. The sound function is now disabled.

### 2) MEMORY FUNCTION

- While the appliance is switched off, press the “MEM” button.
- This appliance can memorise 9 sets of measurement values in chronological order (fig. 5)
- The appliance will switch off automatically after 30 seconds of inactivity. Alternatively, it is possible to press the “ⓘ” button to switch the appliance off manually.

Fig. 5
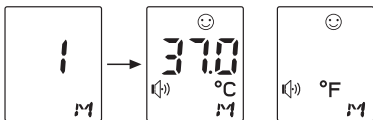
### 3) DISPLAY BACKLIGHTING (for body temperature measurements)

- When the measured temperature is  $<34^{\circ}\text{C}$ , a “LO” message will appear on the display, which will be backlit in red.
- When the temperature measured is  $34^{\circ}\text{C} - 37.1^{\circ}\text{C}$ , the display will be backlit in green.
- When the temperature measured is  $37.2^{\circ}\text{C} - 38^{\circ}\text{C}$ , the display will be backlit in orange.
- When the temperature measured is  $38.2^{\circ}\text{C} - 43.0^{\circ}\text{C}$ , the display will be backlit in red.
- When the temperature measured is  $>43.0^{\circ}\text{C}$ , the display will be backlit red and a “HI” message will appear on the display.

## INSTRUCTIONS FOR USE

### 1) BODY TEMPERATURE MEASUREMENTS

The product must be used to measure human body temperature only.

The temperature of the human body is almost constant, but it can change slightly. It changes throughout the day even in normal conditions of health.

For example, on average:

#### **At night**

Body temperature is lower due to the low level of activity.

#### **In the morning**

When we get up, because the temperature of the room is lower than the temperature in bed, the body's muscles contract and produce heat, which raises the temperature.

#### **At midday**

The human body reaches its highest temperature after lunch.

#### **Mid-afternoon**

Due to physical exertion, the blood glucose level drops and the body's average temperature drops.


#### **In the evening**

As the sun goes down, the temperature of the body drops once again.

This appliance can be used to measure human body temperature on the forehead.

**Fig. 6**

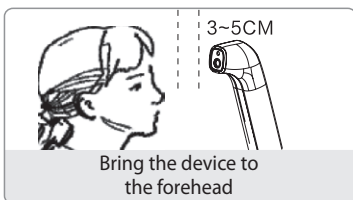
Forehead temperature

## Measuring the temperature at the forehead

- Make sure that there is no perspiration or make-up on the forehead and temples and remember that the forehead measurement can be affected by the use of vasoconstrictor medication and skin irritation.
- Press the "  " key to switch the thermometer on. The LCD display will show the last temperature recorded.
- Direct the thermometer towards the forehead at a distance of 3 - 5 cm
- Press the "  " button for one second. Once the measurement has been taken, a short tone will sound and the value measured will appear on the display (Fig. 7).

**Fig. 7**

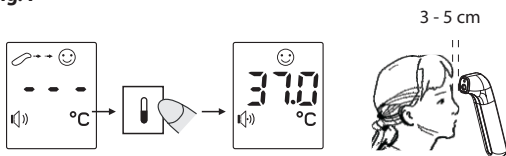When LCD flashing, we are waiting for starting measurement


Measurement result

*Note: If you do not hear the tone, the measurement has not yet been completed. Do not move the thermometer from the forehead yet. (No tone will be heard if the sound function has been disabled).*

## 2) MEASUREMENT OF THE TEMPERATURE OF LIQUIDS/ OBJECTS AND THE ENVIRONMENT

- Hold down the “**I**” button to obtain the object measurement mode.
- When this mode is used, the range of measurement is different to that for body temperature measurements. In this case, the range is 0 °C - 93.2 °C.
- Direct the thermometer at the object, liquid or environment to be measured.
- Press the “**U**” button for one second. Once the measurement has been taken, a short tone will sound and the value measured will appear on the display (Fig. 8).

**Fig. 8**

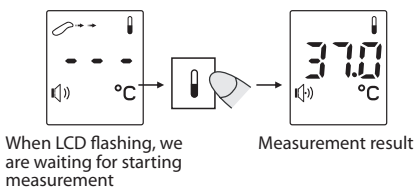*Note: If the tone does not sound, the temperature measurement has not yet been completed. Do not move the thermometer away from the object, liquid or environment yet. (No tone will be heard if the sound function has been disabled).*

## BATTERY REPLACEMENT

When the display shows the “LO” message and the battery symbol is fixed, replace the batteries with batteries of the same type.

1. Open the battery compartment as indicated in figure 3.
2. Remove the exhausted batteries and dispose of them as indicated in the “Disposal” chapter.
3. Insert batteries of the same type and with the same service life, making sure they are positioned as indicated.
4. Dispose of the exhausted batteries in compliance with applicable local regulations.
5. Close the battery compartment.


## MAINTENANCE AND UPKEEP

In order to guarantee accurate measurements, the device should be cleaned after each use.

Use a cotton pad wetted with a 75% rubbing alcohol solution to remove any residual dirt. Take at least 5 seconds to clean the probe. Then clean the other parts of the appliance using a soft dry cloth for at least 10 seconds. After cleaning, make sure there are no visible dirt or residues on the thermometer, and wait for 10 minutes before taking the measurement.

To avoid damaging the instrument or its parts, use soft materials for cleaning.

**Fig. 9**
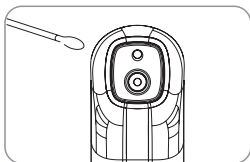
### **Why should the device be sanitised after each measurement?**

This is necessary because infrared thermometers use a highly sensitive technique to measure temperature. Not only is the accuracy of the measurement affected by the presence of dirt and/or dust, but an inadequately sanitised device may also spread infection.

## DISPOSAL

The product packaging is made of recyclable materials. It must be disposed of in compliance with local regulations.

 Dispose of exhausted and completely discharged batteries in the appropriate collection points, at waste collection site for toxic waste or at electronics stores. Battery disposal is a legal requirement.

 For environmental reasons, the appliance must not be discarded with normal household waste at the end of its service life. It must be taken to an appropriate disposal facility. Dispose of the appliance in accordance with European Directive 2012/19/EU regarding waste electrical and electronic equipment. For further information or clarification, contact the local authority responsible for waste disposal.

## AFTER-SALES SERVICE AND WARRANTY

**This appliance is intended for domestic use only and is not suitable for professional purposes.**

- The appliance is covered by a two-year warranty from the time of purchase, which must be proven by proof of purchase. This period complies with applicable regulations ("Consumer Code" Italian Legislative Decree no. 206 of 06/09/2005) and only applies if the consumer is a private subject.
- The warranty covers production defects only and is not valid for damage caused by accidents, incorrect use, negligence or improper use of the product. It does not cover materials subject to normal wear and tear.
- The appliance must not be opened under any circumstance; the warranty is not valid if the product has been opened or tampered with.
- The warranty expires 2 years after purchase; after this time, any technical assistance will be provided on a paid basis. For information on technical assistance under warranty or provided on a paid basis, contact the manufacturer using the details on the last page of this manual. Consumers will not receive any reimbursement for product repairs or replacements performed under warranty. In the event of a malfunction, contact the manufacturer; Do NOT send the appliance directly to Farmamed, **send an e-mail to [service@farmamed.it](mailto:service@farmamed.it)** first. Work performed under warranty (including replacement of all or part of the product) will not extend the duration of the original warranty of the replaced product.
- The manufacturer declines all responsibility for any damage/injury that may be directly or indirectly caused to people, property or domestic animals as a consequence of failure to comply with the indications provided in this instruction manual, in particular the warnings regarding the installation, use and maintenance of the appliance.
- Given its on-going commitment to improving its products, Farmamed reserves the right to modify all or parts of its products in order to meet production requirements, without giving advanced notice, without this entailing any responsibility for Farmamed or its dealers.

## TROUBLE-SHOOTING

In the event of problems during use of the appliance, consult the information below.

Problem	Possible causes	Solutions
The display does not light up when the thermometer is switched on.	The batteries are exhausted.	Replace with new batteries.
	The batteries are inserted the wrong way around.	Insert the batteries as shown on the battery compartment.
The measured temperature is low.	The position used to take the measurement is incorrect.	Measure the temperature correctly as described in the instructions.
	The sensor is dirty.	Remove all dirt before taking the measurement.
Significant temperature changes with continuous measurements.	The measurement interval is too short.	The interval for each measurement should be longer than 10 seconds.
Hi	When the temperature is outside the measurement range, i.e. when the temperature is over 43.0 °C, the LCD display will show the "Hi" message.	The device cannot be used.
Lo	When the temperature of the target object is outside the measurement range, i.e. when the temperature is under 34 °C, the LCD display will show the "Lo" message.	The device cannot be used.
Er.H	This product's maximum operating temperature is 43 °C. When the temperature exceeds this value, the "Er.H" error message will appear on the LCD display.	The device cannot be used.
Er.L	This product's lowest operating temperature is 10°C. When the temperature drops below this value, the "Er.L" error message will appear on the LCD display.	The device cannot be used.
Err	When the ambient temperature changes quickly, the "Err" message will appear during measurements and the device will switch off automatically.	When the "Err" message appears, put the product down and leave it for 30 minutes before taking another measurements.

## PRODUCT SPECIFICATION

**Product name:**

Infrared non-contact Thermometer

**Product dimensions:** 34 x 160 x 50 mm

**Weight:** about 68g (batteries excluded)

**Measurement range (body temperature):** 34.0°C-43.0°C.  
(93.2°F-109.4°F)

**Measurement range (liquids, objects and environments):**  
0°C-93.2°C (32°F-199.7°F)

**Precision:** (35.0°C-42.0°C) ± 0.2°C, (95.0°F-107.9°F) ± 0.4°F, ultra  
temperatura ± 0.3°C.

**Working conditions:**

- **Temperature:** 10.0°C-40.0°C(50.0°F-104.0°F)
- **Maximum relative humidity:** 15%RFI-93%RFI
- **Atmospheric pressure:** 70kPa-106kPa

**Conditions for storage and transportation:**

- **Temperature:** -25°C-55°C (-13°F-131°F),
- **Maximum relative humidity:** 0%RFI-93%RFI
- **Atmospheric pressure:** 50.0-106.0kPa

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Shenzhen Pango Medical Electronics Co.,Ltd
Lotus NL B.V.
**Farmamed S.r.l.**