



**WALMAN**  
INSTRUMENTS

**BQE**  
BACK QUALITY ERGONOMICS  
STOOLS & CHAIRS  
DUTCH DESIGN<sup>®</sup>




## WM ESB

The ESB is developed by a Dutch PT. The saddle/bicycle shaped seat allows for higher positioning and facilitates a better posture. The convex backrest fits snugly in the concave curvature of the lower back and prevents flattening of the lumbar area. The mold injected foam has a 3d shape preventing circulation problems on the legs and allowing for almost full leg extension.



## WM CST

The CST is developed by a Dutch PT. The triangular shaped seat allows for higher positioning and facilitates a better posture. The convex backrest fits snugly in the concave curvature of the lower back and prevents flattening of the lumbar area. The mold injected foam has a 3d shape preventing circulation problems on the legs and allowing for almost full leg extension.

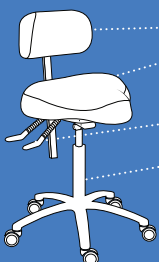

**FOAM & UPHOLSTREY:** The mold injected foam is the best quality foam. The vinyl colo is standard in black has special layer of and makes it easy to clean and durable. Additional colors are available at additional cost.

**MECHANISM:** The mecanism offers a seat tilt, backrest adjustment and height adjustment for the stool.

**BASE:** The five-legged base ensures maximum stability, the 60mm casters creates easy movement around the patient (braked casters or gliders are available at additional cost).

**DIMENSIONS:** Base diameter 23" – Seat 17" x 16" – Back 11" x 7" – Seat height 20"-28" – Weight 33lbs.