
# FX-920



# FX-920

# The first full-power tilt chair.

With the Model FX-920 Tilt Chair, the patient maintains his position as the chair is reclined. The practitioner may operate all functions while maintaining contact with patient and remaining seated.

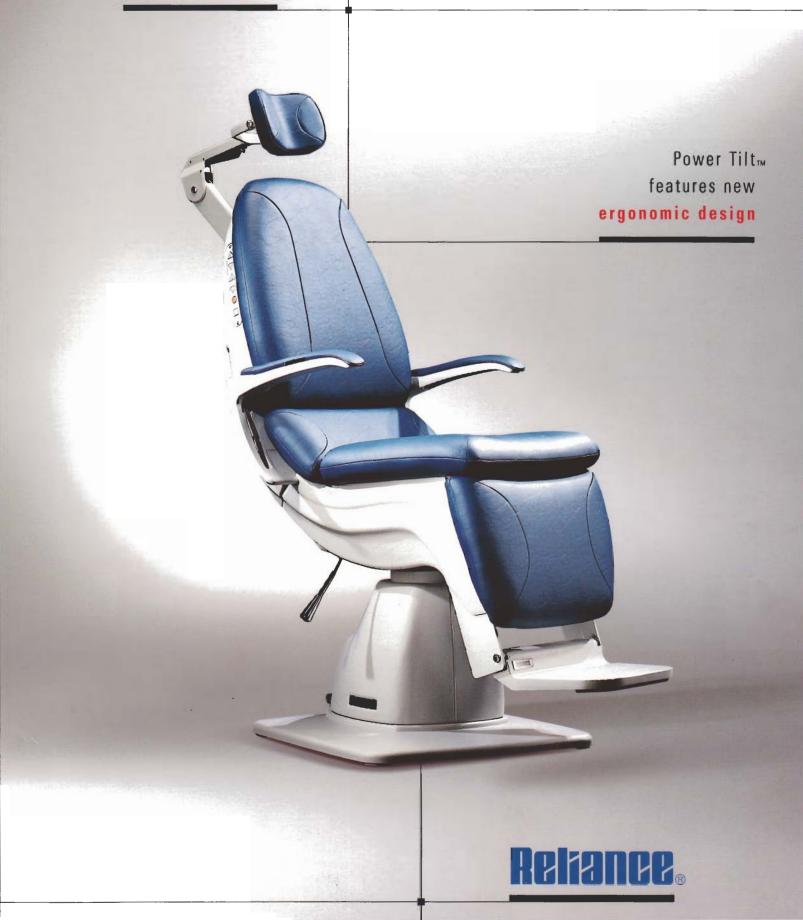

# Power Reclining Top

Use membrane finger switches and corded footswitch to raise, lower or tilt chair into any position without breaking contact with patient. Memory programming eliminates wasted movements. An auto-return button automatically returns the patient to the low, upright position. Patient maintains position as chair reclines (tilts).

# **Dual Side Operation**

For comfort and efficiency, chair design enables practitioner to execute any chair function from either side and still maintain contact with patient.

## Power Hydraulic Base

Backlit membrane switches control the raise, lower and auto down functions of the FX-920 examination chair. Membrane switches last 10-times longer than their mechanical counterparts, and they increase safety by eliminating high voltages.

#### **Dual Rotation Locks**

Easily accessible levers located on both sides of seat enable practitioner to swivel chair 350° without changing position.

# Large Hinged Footstep

The 8"x 13" footstep can support patient's full weight, and offers safe ingress and egress. Tilts out of way when desired

#### Headrest

Concave upholstered headrest easily adjusts, locks firmly in place and unlocks with a single hand. Wide range of vertical and axial positions. Provides superior head and neck support.

#### Armrests

Upholstered armrests lock in place for patient security and comfort. Release mechanism unlocks them for convenient ingress and egress from either side of chair.

### Finishes and Upholstery

Hand-sewn upholstery with high density foam holds its shape for years. Cushions remove for rapid cleaning and replacement. Match any office or treatment room décor with our selection of 14 easy-to-maintain vinyl colors. Base is available in #20 Pearl. All materials are fire retardant. To request a color selection card, contact your dealer or call Reliance at 1-800-735-0357.

#### **Dimensions and Specifications**

Top and base are shipped together for easy installation. Choice of two base heights: HIGH: 22" lowered, 33 1/2" raised; LOW: 20" lowered, 29" raised. Base: 22 1/8" x 24 1/4". Width: 24 1/4". Depth upright (excluding headrest): 41". Reclined: 54". Seat: 20" (wide) x 15" (deep). Distance between armrests: 19". Shipping weight: 355 lbs. Electrical: (Standard) 110 volts, 60 Hz, 8.0 amperes; (Optional) 230 volts, 50 Hz, 3.0 amperes. UL Approved. Patent pending.

